

WIDEX CLEAR™ 440

STAY IN SYNC WITH LIFE

USER INSTRUCTIONS

**THE CLEAR™ 440 SERIES**

C4-PA model

RIC

Receiver-in-canal

**WIDEX**®

HIGH DEFINITION HEARING

## **YOUR HEARING AID**

*(To be filled out by the hearing care professional)*

**Date:** \_\_\_\_\_

**Battery size:** \_\_\_\_\_

# CONTENTS

<b>SYMBOLS</b> .....	<b>4</b>
<b>THE HEARING AID</b> .....	<b>5</b>
Right/left identification.....	8
Acoustic indicators.....	8
The battery.....	9
Inserting the battery.....	9
Low battery indication.....	11
Turning the hearing aid on and off.....	12
Positioning the hearing aid and ear-tip.....	13
Removing the hearing aid and ear-tip.....	14
Volume adjustment.....	15
Lost partner alarm.....	16
Listening programs.....	17
Switching between the listening programs.....	19
Zen+.....	19
Using a telephone.....	20
Cleaning.....	21
The hearing aid.....	22
<b>ACCESSORIES</b> .....	<b>23</b>
<b>IN CASE OF MALFUNCTION</b> .....	<b>24</b>
<b>CARING FOR YOUR HEARING AID</b> .....	<b>26</b>
<b>WARNINGS</b> .....	<b>27</b>
<b>ADVICE</b> .....	<b>30</b>
<b>REGULATORY INFORMATION</b> .....	<b>31</b>

# SYMBOLS

The following symbols will be used throughout the manual:



Messages with this heading indicate serious adverse reactions, potential safety hazards and inadequate performance of device.



Messages with this heading indicate/include information regarding any special care to be exercised.




 Not for general waste.

# THE HEARING AID

The illustration below shows the hearing aid without the ear-set. The choice of ear-set solution depends on your specific needs. Please refer to the separate ear-set user manual.

1. Microphone openings
2. On/off function
3. Nail grip
4. LED / Right/left identification

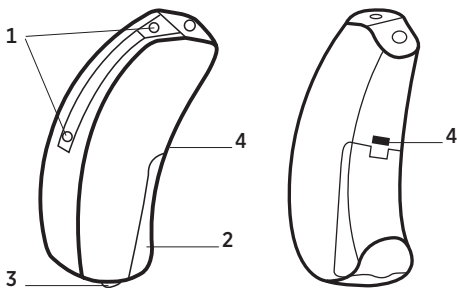

## NOTE

In addition to these user instructions, a separate user manual is provided describing the various ear-set solutions available for your hearing aid: “Ear-sets for Widex BTE hearing aids”.



### **WARNING**

This booklet and the manual “Ear-sets for Widex BTE hearing aids” contain important information and instructions. Read these booklets carefully before you start using the hearing aid.

## NOTE

Your hearing aid, ear-set and accessories may not look exactly as illustrated in this booklet. We also reserve the right to make any changes considered necessary.

## **Intended use**

The CLEAR hearing aids are intended as air conduction amplification devices to be used in everyday listening environments. The CLEAR hearing aids may be provided with the Zen program intended to provide a relaxing sound background (i.e. music/noise source) for adults who desire to listen to such a background in quiet.

## **Indications for use**

The devices are indicated for individuals with a range of hearing loss severity from slight (16 dB HL to 25 dB HL) to moderately severe (85 dB HL) ) and all hearing loss configurations.


They are to be programmed by licensed hearing health-care professionals (audiologists, hearing aid specialists, otolaryngologists) who are trained in hearing (re)habilitation as well as tinnitus management.

### **Description of device**

Your CLEAR™ hearing aid is used with an ear-set which consists of an earwire and an ear-tip, in which the receiver is placed. The hearing aid uses a proprietary wireless technology, WidexLink, to enable communication between the left and the right hearing aids, as well as between the hearing aids and the DEX accessories.

Your hearing aid may be provided with an optional listening program called Zen. It makes musical tones (and sometimes a rushing noise) in the background. These sounds are shaped according to your hearing loss.

## Right/left identification

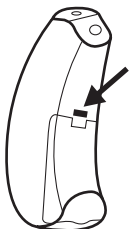

The arrow shows the position of the identification mark (red mark = right and blue mark = left).

## Acoustic indicators

The hearing aid may be set to produce a signal to indicate the use of certain functions. The signal may be a spoken message or tones adjusted by the hearing care professional in consultation with you. The signal can also be deactivated.

Adjusting volume via remote control		<b>Tone*</b>	Off
Confirming RC program key use		<b>Clicking sound</b>	Off
Changing program via remote control	<b>Message</b>	Tones	Off
Starting up the hearing aid	<b>Message</b>	Tone	Off
Warning about low battery	<b>Message</b>	4 tones	Off
Warning about lost partner	Message		<b>Off</b>
Service reminder	Message		<b>Off</b>

\* **Bold** indicates standard setting



## The battery

We recommend **zinc-air batteries**. Use a **size 10** battery for the hearing aid.

To obtain replacement batteries, please consult your hearing care professional. It is important to take note of the expiration date and the recommendations on the battery pack regarding disposal of used batteries. Otherwise there is a risk that the battery life may be reduced.



### Inserting the battery

Before inserting a new battery into the hearing aid, remember to remove the adhesive tab. Once the tab has been removed, the battery will start functioning after a few seconds.




### CAUTION

Do not use batteries if there is a sticky residue from the tab or other unwanted substance, as this can cause the hearing aid to malfunction.



Use the nail grip to gently swing the battery drawer open. Do not press the drawer beyond the open position.

Place the battery in the drawer as illustrated, so that the plus (+) sign on the battery is facing towards you when you hold the hearing aid as shown. You can use the battery magnet provided to steer the battery into place.

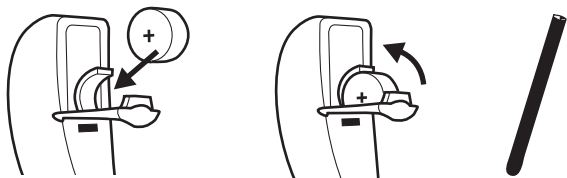

If the battery drawer does not close easily, the battery is incorrectly inserted.

When changing battery, it is a good idea to hold the hearing aid over a table.

Always place the battery in the actual battery drawer as illustrated and not in the open space inside the hearing aid.

## Low battery indication

An acoustic indicator will sound when the battery is nearly exhausted, unless this function has been disabled (see page 8). Exactly how long the hearing aid will function after the low battery indication varies from one case to another. We recommend that you always have a spare battery with you.



### WARNING

**Never** leave an exhausted battery in the hearing aid. Exhausted batteries may leak, damaging the hearing aid.




### WARNING


Your hearing aid may stop functioning, for instance if the battery is exhausted. You should be aware of this possibility, in particular when you are in traffic or are otherwise dependent on warning signals.

## Turning the hearing aid on and off

The battery drawer of the hearing aid also functions as the on/off switch.

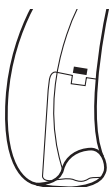

Close the battery drawer to turn on the hearing aid. An acoustic indicator will indicate that the hearing aid has been switched on unless this function has been disabled.

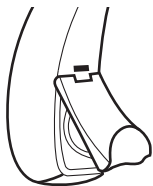



Open the battery drawer to turn off the hearing aid.

Please remember to turn off the hearing aid when it is not in use. Remove the battery if the hearing aid will not be used for several days.

## Positioning the hearing aid and ear-tip

Place the hearing aid behind the ear, so that the earwire rests comfortably on the ear, close to your head.

Insert the ear-tip into the ear canal while holding the lower part of the earwire.

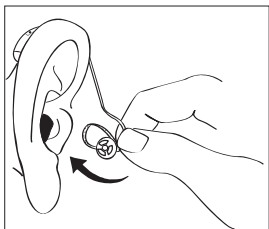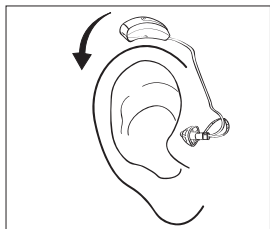


The illustrations show an open ear-tip. For further information on ear-tip types, anchors and procedures, see the separate ear-set user manual provided.

## Removing the hearing aid and ear-tip

Remove the hearing aid from its position behind the ear.

Carefully pull the ear-tip out of the ear canal, while holding the lower part of the earwire. If the ear-tip is provided with an extraction cord, take hold of this and carefully pull the ear-tip out of the ear canal.

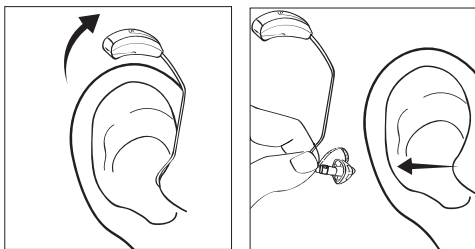

## Volume adjustment

The hearing aid volume is adjusted automatically in accordance with your sound environment.

If you have a remote control, you can also adjust the volume manually or mute your hearing aid. Any adjustment of the volume setting will be cancelled when your hearing aid is turned off, or when you change program.

In the default setting, any volume adjustment via the remote control will affect both hearing aids.

To alter the increments in which the volume is increased or decreased, consult your hearing care professional.

If you have a remote control, please also refer to the user instructions for this.

If you would like any further information, consult your hearing care professional.



### **CAUTION**

If the volume in the hearing aid is generally too loud or too weak, or the reproduced sounds are distorted, consult your hearing care professional.

### **To mute the hearing aid with a remote control:**

Keep pressing the volume down key on the remote control after the long beep-tone has sounded and until it stops. Pressing one of the volume keys briefly will bring back the sound.

### **Lost partner alarm**

**(Note: Only applicable if this feature has been activated by your hearing care professional).**

If one of your hearing aids falls off your ear, and/or the battery is exhausted, a tone signal and a spoken message will sound two times in the other hearing aid, and a diode will flash on both hearing aids.

When the lost hearing aid is repositioned, and/or the battery is exchanged, the diodes stop flashing.



## Listening programs

By default, your hearing aid has one listening program. If you have a remote control, your hearing aid may be provided with five of the following programs, plus the special Zen program called Zen+.

Please also refer to the user instructions for your remote control.

<input type="checkbox"/> <b>Master</b>	Standard
<input type="checkbox"/> <b>Music</b>	For listening to music
<input type="checkbox"/> <b>TV</b>	For listening to the TV
<input type="checkbox"/> <b>Comfort</b>	Attenuates background noise
<input type="checkbox"/> <b>Reverse focus</b>	Focus on sounds from behind
<input type="checkbox"/> <b>Audibility Extender</b>	Makes high-frequency sounds audible
<input type="checkbox"/> <b>Zen</b>	Generates different kinds of tones
<input type="checkbox"/> <b>Phone</b>	For listening to the telephone
<b>Compound programs</b>	The Master program in one ear and Zen, Audibility Extender, Reverse focus or Phone in the other
<input type="checkbox"/> _____	
<input type="checkbox"/> _____	
<input type="checkbox"/> _____	
<input type="checkbox"/> _____	
<input type="checkbox"/> <b>Zen+</b>	Special Zen program with up to three different Zen styles

If your needs change over time, the combination of listening programs can be altered.

The Zen program may be used with amplification so that both the surrounding sounds and the Zen tones are heard together. Or, it may be used alone (without amplification) in quiet when you are not required to hear surrounding sounds.

The duration of the Zen play time can be set by your hearing care professional according to your needs.

### **Benefits**

The Zen program may provide a relaxing listening background for some people. When the Zen program is used in a tinnitus management program, its wearer may experience some relief from tinnitus.



#### **CAUTION**

Use of the different Zen programs may interfere with hearing surrounding sounds including speech. The programs should not be used when hearing such sounds is important. Switch the hearing aid to a non-Zen program in those situations.



#### **CAUTION**

If you perceive decrease in loudness, tolerance of sounds, speech not as clear, or worsening tinnitus, contact your hearing care professional.

## Switching between the listening programs

Each time you change program via the remote control, an acoustic indicator will sound unless this function has been disabled.

**Program 1:** Message or one short beep

**Program 2:** Message or two short beeps

**Program 3:** Message or three short beeps

**Program 4:** Message or one long and one short beep

**Program 5:** Message or one long and two short beeps


**Zen+:** Tone or message

In the default setting, both your hearing aids will be affected when you change programs.

### Zen+

This program is accessed via a long press (two secs.) on the RC-DEX program key. A short press then allows you to cycle through the available Zen styles. You can exit Zen+ by pressing and holding the program key for two seconds.

## Using a telephone

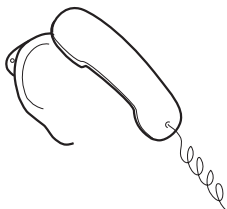


When using a telephone, we recommend that you hold the telephone against your head at an angle above your ear, rather than directly against the ear. If the sound is not optimal, try moving the telephone earpiece to a slightly different position.

## Cleaning


The following cleaning accessories are available for the hearing aid and ear-set\*. For cleaning the ear-set, see the user manual “Ear-sets for Widex BTE hearing aids.”

1. Cloth
2. Long wax removing tool
3. Short wax removing tool
4. Brush

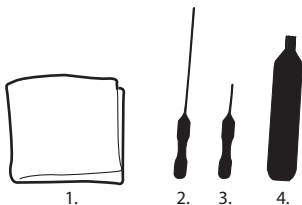

Contact your hearing care professional if you need additional supplies of cleaning accessories.

*\* Selection depends on ear-set type.*

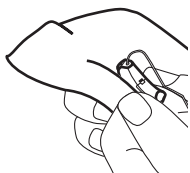


### The hearing aid

Clean the hearing aid after use with the cloth provided.



#### **WARNING**

**Never** use water or cleaning solutions to clean the hearing aid, as this may cause it to malfunction.

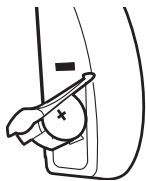

When the hearing aid is not in use, the battery drawer should be left open, so as to ventilate the hearing aid and allow it to dry.

## ACCESSORIES

A variety of assistive listening devices are available for your hearing aid.

- RC-DEX
- TV-DEX
- M-DEX

For help in determining whether you could benefit from a remote control or other accessories, please consult your hearing care professional.

## IN CASE OF MALFUNCTION

The following pages contain some quick advice on what to do if your hearing aid stops working or performs unsatisfactorily. If the problems persist, contact your hearing care professional for assistance.

See the “Ear-sets for Widex BTE hearing aids” user manual for information specific to your ear-set.

<b>Problem</b>	<b>Potential cause</b>	<b>Solution</b>
The hearing aid is completely silent	The hearing aid is not turned on	Make sure the battery drawer is closed
	The battery does not work	Insert a new battery in the hearing aid
The hearing aid volume is not powerful enough	The battery is nearly exhausted	Insert a new battery in the hearing aid
	Your ear is blocked by earwax	Contact your ENT doctor/physician
	Your hearing may have changed	Contact your hearing care professional
The hearing aid whistles continuously	Your ear is blocked by earwax	Contact your ENT doctor/physician
Your two hearing aids are not working in synchrony	The connection between the hearing aids is lost	Turn the hearing aids off and on again



<b>Problem</b>	<b>Potential cause</b>	<b>Solution</b>
The hearing aids do not respond with a corresponding change in volume or program to the DEX	<ul style="list-style-type: none"> <li>a. The DEX is used beyond the transmission range</li> <li>b. Strong electromagnetic interference in the vicinity</li> <li>c. The DEX and the CLEAR hearing aids are not matched</li> </ul>	<ul style="list-style-type: none"> <li>a. Move the DEX closer to the CLEAR hearing aids</li> <li>b. Move away from known source of EM interference</li> <li>c. Check with hearing care professional to make sure DEX is matched with CLEAR hearing aids</li> </ul>
You hear “interrupted” speech (on and off) from the hearing aids or no speech (muted) from the transmitting ear.	<ul style="list-style-type: none"> <li>a. The battery in one of the hearing aids has expired</li> <li>b. Strong electromagnetic interference in the vicinity</li> </ul>	<ul style="list-style-type: none"> <li>a. Replace battery in one or both hearing aids</li> <li>b. Move away from known sources of interference</li> </ul>

## CARING FOR YOUR HEARING AID

The hearing aid is a valuable object and should be treated with care. Here are some things you can do to prolong the life of your hearing aid:



### CAUTION

- Turn off your hearing aid when it is not in use. Remove the battery if the hearing aid will not be used for several days.
- When the hearing aid is not in use, keep it in its case in a cool, dry location out of reach of children and pets.
- Do not expose the hearing aid to extreme temperatures or high humidity. Make sure to dry the hearing aid thoroughly after heavy perspiration such as that which may occur during intense physical activity, e.g. playing sports.
- Avoid dropping your hearing aid – perform cleaning and battery changes while holding the hearing aid above a soft surface.
- Do not wear your hearing aid while in the shower or swimming, or when using a hair dryer, perfume, hair and body sprays or gels such as suntanning lotions or creams.
- In environments with high humidity, a drying kit may be used daily to reduce the amount of moisture inside the hearing aid. See the instructions enclosed with the drying kit.

## WARNINGS



### WARNING

**Hearing aids and batteries can be dangerous if swallowed or used improperly. Swallowing or improper use can result in severe injury or even fatalities. In case of ingestion, contact your physician immediately.**

- Keep hearing aids and their parts, accessories and batteries out of reach of children and anyone else who might swallow such items or otherwise cause injury to themselves. Do not change batteries in front of them and do not let them see where you keep your battery supply. Discard used batteries carefully.
- Batteries are very small and can easily be mistaken for pills or the like. Never put a battery or hearing aid in your mouth for any reason as you may risk swallowing it.
- Risk of explosion if battery is replaced by an incorrect type or recharged. Dispose of used batteries according to the instructions.
- Never allow others to wear your hearing aid, as this could cause permanent damage to their hearing.
- When selecting a listening program, please remember that there are situations in which it is particularly important to be able to hear the surrounding sounds (e.g. traffic, warning signals).
- The hearing aid is made of modern non-allergenic materials. Nonetheless, in rare cases skin irritation may occur. If you notice skin irritation in or around your ear or ear canal, contact your hearing care professional.



## WARNING

- Be aware that the use of any type of hearing aid may involve a slightly increased risk of infection in the ear canal. An infection can arise as a result of inadequate ventilation of the ear. We therefore recommend that you remove the hearing aid and ear-set from your ear at night to allow the ear canal to be ventilated. Make sure that you clean and inspect your hearing aid and ear-set as required. If an infection occurs, you should seek medical attention and contact your hearing care professional for advice on how to disinfect the various hearing aid parts. Do not under any circumstances use alcohol, chlorine or similar substances.
- Do not use Widex hearing aids in mines or other areas with explosive gases.
- **Do not wear your hearing aid during radiation, X-rays, MRIs, CT or other medical treatments and scans.** The emissions from these procedures as well as from other types of radiation, such as that in a microwave oven, can damage your hearing aid. Radiation from, for example, room surveillance equipment, burglar alarms and mobile phones is weaker and will not damage the hearing aid, but may create audible interference.



### **CAUTION**

- Your hearing aid has been tested for interference according to international standards. Nevertheless, it is possible that unforeseen interference may occur in the hearing aid due to electromagnetic radiation from other products such as alarm systems, room surveillance equipment and mobile phones.
- Although your hearing aid has been designed to comply with the most stringent international electromagnetic compatibility standards, the possibility cannot be excluded that it may cause interference with other equipment, such as medical devices.
- Never try to open or repair the hearing aid yourself.

## ADVICE

### NOTE

- The hearing aid will not restore normal hearing and will not prevent or improve a hearing loss resulting from organic conditions. However, the hearing aid can help you to make the best possible use of your remaining hearing ability. You should also bear in mind that it can take time to get used to a new hearing aid and new sounds.
- In most cases, using the hearing aid infrequently will not permit you to gain full benefit from it.
- The use of a hearing aid is only part of hearing habilitation and may need to be supplemented by auditory training and instruction in lipreading.
- The use of hearing aids increases the risk of accumulation of earwax. Contact your physician/ENT doctor if you suspect that a plug of earwax has accumulated in your ear. Earwax may not only reduce your own hearing but also the effect of the hearing aid considerably. It is a good idea to ask your physician to clean your ears a couple of times a year.

# REGULATORY INFORMATION

FCC ID: TTY-C4PA

IC: 5676B-C4PA

## Federal Communications Commission Statement

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

### NOTE:

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

**NOTE:**

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

Changes or modifications to the equipment not expressly approved by Widex could void the user's authority to operate the equipment.

## **Industry Canada Statement / Déclaration d'industrie Canada**

Under Industry Canada regulations, this radio transmitter may only operate using an antenna of a type and maximum (or lesser) gain approved for the transmitter by Industry Canada.

To reduce potential radio interference to other users, the antenna type and its gain should be so chosen that the equivalent isotropically radiated power (e.i.r.p.) is not more than that necessary for successful communication.

This device complies with Industry Canada licence-exempt RSS standard(s). Operation is subject to the following two conditions:

- (1) this device may not cause interference, and
- (2) this device must accept any interference, including interference that may cause undesired operation of the device.



Conformément à la réglementation d'Industrie Canada, le présent émetteur radio peut fonctionner avec une antenne d'un type et d'un gain maximal (ou inférieur) approuvé pour l'émetteur par Industrie Canada.

Dans le but de réduire les risques de brouillage radioélectrique à l'intention des autres utilisateurs, il faut choisir le type d'antenne et son gain de sorte que la puissance isotrope rayonnée équivalente (p.i.r.e.) ne dépasse pas l'intensité nécessaire à l'établissement d'une communication satisfaisante.

Le présent appareil est conforme aux CNR d'Industrie Canada applicables aux appareils radio exempts de licence.

L'exploitation est autorisée aux deux conditions suivantes :

- (1) l'appareil ne doit pas produire de brouillage, et
- (2) l'utilisateur de l'appareil doit accepter tout brouillage radioélectrique subi, même si le brouillage est susceptible d'en compromettre le fonctionnement.

CE 0459

Hereby, Widex A/S declares that this C4-PA is in compliance with the essential requirements and other relevant provisions of Directive 1999/5/EC.

A copy of the Declaration of Conformity can be found at: <http://www.widex.com>



Hearing aids, accessories and batteries should not be disposed of with normal household waste. Please consult your national Widex distributor for advice on how to dispose of these items.



**CE** 0459

**Manufacturer**

**WIDEX A/S** • Nymoellevej 6, DK-3540 Lyngø  
Denmark • [www.widex.com](http://www.widex.com)



Printed by HTO / 2011-07  
9 514 0134 001 #02